


HXP50-MECA Hexapod

6-Axis-Parallel Kinematic Positioning System

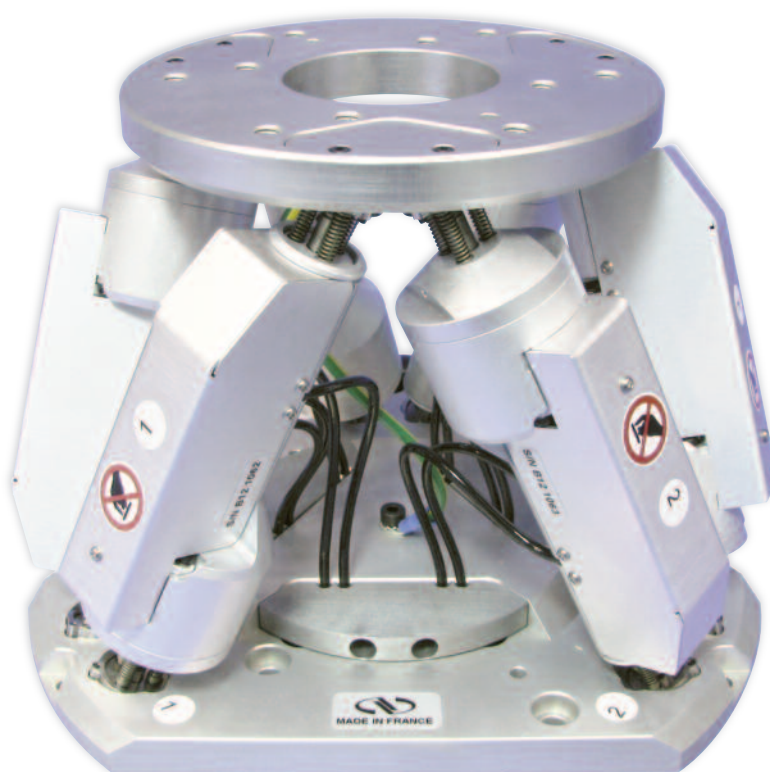

Newport®

Experience | Solutions

USER'S MANUAL

Precision Motion – **Guaranteed™**

Warranty

Newport Corporation warrants this product to be free from defects in material and workmanship for a period of 1 year from the date of shipment. If found to be defective during the warranty period, the product will either be repaired or replaced at Newport's discretion.

To exercise this warranty, write or call your local Newport representative, or contact Newport headquarters in Irvine, California. You will be given prompt assistance and return instructions. Send the instrument, transportation prepaid, to the indicated service facility. Repairs will be made and the instrument returned, transportation prepaid. Repaired products are warranted for the balance of the original warranty period, or at least 90 days.

Limitation of Warranty

This warranty does not apply to defects resulting from modification or misuse of any product or part.



CAUTION

Please return equipment in the original (or equivalent) packing.

You will be responsible for damage incurred from inadequate packaging if the original packaging is not used.

CAUTION

Warranty does not apply to damages resulting from:

- **Incorrect usage:**
 - Load on the Hexapod greater than maximum specified load.
 - Carriage speed higher than specified speed.
 - Improper grounding.
 - Connectors must be properly secured.
 - When the load on the Hexapod represents an electrical risk, it must be connected to ground.
 - Excessive or improper cantilever loads.
- **Modification of the Hexapod or any part thereof.**

This warranty is in lieu of all other warranties, expressed or implied, including any implied warranty of merchantability or fitness for a particular use. Newport Corporation shall not be liable for any indirect, special, or consequential damages.

No part of this manual may be reproduced or copied without the prior written approval of Newport Corporation.

This manual has been provided for information only and product specifications are subject to change without notice. Any changes will be reflected in future printings.

Table of Contents

| | |
|---|--------|
| Warranty | ii |
| EC Declaration of Conformity..... | v |
| Definitions and Symbols..... | vi |
| Warnings | vii |
| Cautions | viii |
| 1.0 — Introduction | 1 |
| 1.1 Remove the Locking System..... | 2 |
| 2.0 — Description | 3 |
| 3.0 — Characteristics | 4 |
| 3.1 Mechanical Specifications | 4 |
| 3.2 Max. Cantilever Distance of the Load..... | 4 |
| Horizontal Base Plate | 4 |
| Base Plate Upside-Down..... | 4 |
| Vertical Base Plate | 5 |
| Base Plate at Any Position | 5 |
| Horizontal Base Plate - Lateral Force | 5 |
| 3.3 Weight..... | 5 |
| 4.0 — Connection to the HXP50-ELEC Controllers | 6 |
| 4.1 Warnings on Controllers | 6 |
| 4.2 Connection..... | 7 |
| 4.3 Cables | 7 |
| 5.0 — Dimensions | 9 |
| 6.0 — Maintenance | 10 |
| 6.1 Maintenance | 10 |
| 6.2 Repair | 10 |
| Service Form | 11 |

EC Declaration of Conformity

HXP50-MECA Series



EC Declaration of Conformity

following Annex II-1A
of Directive 2006/42/EC on machinery

The manufacturer:

MICRO-CONTROLE Spectra-Physics,
9, rue du Bois Sauvage
91055 Evry CEDEX, FRANCE

Hereby declares that the machinery:

- Description: " HXP50-MECA "
- Function: 6-Axis Kinematic Positioning System
- Models: HXP50-MECA ; HXP50-MECAV6

– the technical file of which was compiled by:

Mr Dominique DEVIDAL, Quality Director,
MICRO-CONTROLE Spectra-Physics, Zone Industrielle - B.P.29
F-45340 Beaune La Rolande France

– complies with all the relevant provisions of the Directive 2006/42/EC on machinery.
– complies with all the relevant provisions of the Directive 2014/30/EU relating to electro-magnetic compatibility.

– was designed and built in accordance with the following harmonised standards:

- NF EN 61326-1:2013 « Electrical equipment for measurement, control and laboratory use – EMC requirements – Part 1: General requirements »
- NF EN 55011:2010/A1:2011 Class A

– was designed and built in accordance with the following other standards:

- NF EN 61000-4-2
- NF EN 61000-4-3
- NF EN 61000-4-4
- NF EN 61000-4-5
- NF EN 61000-4-6

ORIGINAL DECLARATION

Done in Beaune La Rolande on 26 June 2015
Dominique DEVIDAL
Quality Director

DC1-EN rev:A

Definitions and Symbols

The following terms and symbols are used in this documentation and also appear on the product where safety-related issues occur.

General Warning or Caution



The exclamation symbol may appear in warning and caution tables in this document. This symbol designates an area where personal injury or damage to the equipment is possible.

The following are definitions of the Warnings, Cautions and Notes that may be used in this manual to call attention to important information regarding personal safety, safety and preservation of the equipment, or important tips.



WARNING

Warning indicates a potentially dangerous situation which can result in bodily harm or death.



CAUTION

Caution indicates a potentially hazardous situation which can result in damage to product or equipment.

NOTE

Note indicates additional information that must be considered by the user or operator.

European Union CE Mark



The presence of the CE Mark on Newport Corporation equipment means that it has been designed, tested and certified as complying with all applicable European Union (CE) regulations and recommendations.

Warnings and Cautions



ATTENTION

This stage is a Class A device. In a residential environment, this device can cause electromagnetic interference. In this case, suitable measures must be taken by the user.

Warnings



WARNING

The motion of objects of all types carries potential risks for operators. Ensure the protection of operators by prohibiting access to the dangerous area and by informing the personnel of the potential risks involved.

WARNING

Do not use this Hexapod when struts are emitting smoke or is unusually hot to the touch or are emitting any unusual odor or noise or is in any other abnormal state.

Stop using the Hexapod immediately, switch off struts power and then disconnect the HXP controller from power supply.

After checking that smoke is no longer being emitted contact your Newport service facility and request repairs. Never attempt to repair the Hexapod yourself as this can be dangerous.

WARNING

Make sure that this Hexapod is not exposed to moisture and that liquid does not get into the Hexapod.

Nevertheless, if any liquid has entered the Hexapod, switch off struts power and then disconnect the HXP controller from power supply.

Contact your Newport service facility and request repairs.

WARNING

Do not insert or drop objects into this Hexapod, this may cause an electric shock, or lock the drive.

Do not use this Hexapod if any foreign objects have entered the Hexapod. Switch off struts power and then disconnect the HXP controller from power supply.

Contact your Newport service facility for repairs.

WARNING

Do not place this Hexapod in unstable locations such as on a wobbly table or sloping surface, where it may fall or tip over and cause injury.

If this Hexapod has been dropped or the case has been damaged, switch off struts power and then disconnect the HXP controller from power supply.

Contact your Newport service facility and request repairs.

WARNING

Do not attempt to modify this Hexapod; this may cause an electric shock or downgrade its performance.

WARNING

Do not exceed the usable depth indicated on the mounting holes (see section “Dimensions”). Longer screws can damage the mechanics or cause a short-circuit.

Cautions

CAUTION

Do not place this Hexapod in a hostile environment such as X-Rays, hard UV,... or in any vacuum environment.

CAUTION

Do not use the top plate for handling. This may cause it to malfunction.

CAUTION

Do not manipulate struts. This may cause it to malfunction.

CAUTION

Do not place this Hexapod in a location affected by dust, oil fumes, steam or high humidity. This may cause an electric shock.

CAUTION

Do not leave this Hexapod in places subject to extremely high temperatures or low temperatures. This may cause an electric shock.

- Operating temperature: +10 to +35 °C
 - Storage temperature: -10 to +40 °C (in its original packaging)
-

CAUTION

Do not move this Hexapod if struts power is on.

Make sure that the cable to the electronics is disconnected before moving the Hexapod. Failure to do so may damage the cable and cause an electrical shock.

CAUTION

Be careful that the Hexapod is not bumped when it is being carried. This may cause it to malfunction.

CAUTION

When handling this Hexapod, always unplug the equipment from the power source for safety.

CAUTION


Contact your Newport service facility to request cleaning and specification control every year.

6-Axis-Parallel Kinematic Positioning System

HXP50-MECA Hexapod

1.0 Introduction

This manual provides operating instructions for the Hexapod that you have purchased.

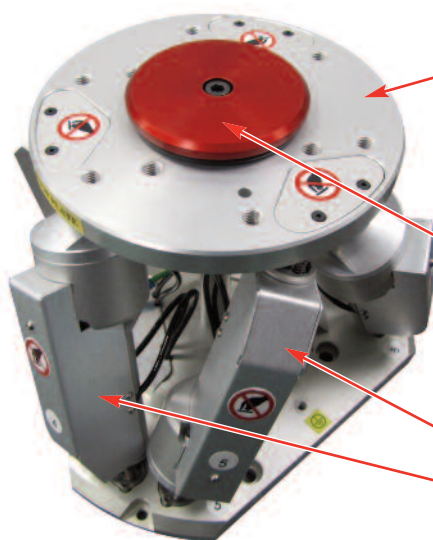**CAUTION**

Do not use the top plate or struts for handling.

**CAUTION**

Remove the locking system before using.

**CAUTION**

Do not manually turn struts.

**CAUTION**


Please read all chapters of this manual carefully before operating your Hexapod.

**NOTE**

This product complies with the RoHS directive (Restriction of Hazardous Substances).

1.1 Remove the Locking System

Hexapod ball studs can be damaged if they are subject to pulling. This can happen while handling or transportation of the HXP50-MECA.

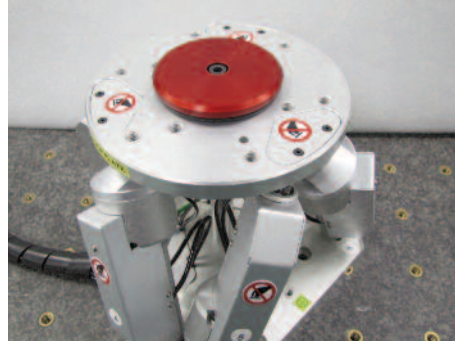

To avoid damages, the HXP50-MECA Hexapod is delivered with a locking system that includes a center rod and a clamping wheel.



CAUTION

Remove the locking system before operating the HXP50-MECA Hexapod.

- 1 Remove the clamping wheel.



- 2 Remove the center rod.



2.0 Description

The HXP50-MECA hexapod is a parallel kinematic motion device that provides six degrees of freedom: X, Y, Z, pitch, roll, and yaw. Hexapods are well known as effective solutions for complex motion applications that demand high load capacity and accuracy in up to six independent axes. Yet, they still have the image of being complex and highly priced. With the HXP50, Newport continues to introduce Hexapods that are not only affordable but also extremely easy to use.


The HXP50-MECA is driven by six DC servo motor driven actuators with encoder feedback at the leadscrew nut, providing precise MIM, low backlash and fast speed. The quality of the actuators has obviously a strong impact on the overall motion performance of the HXP50-MECA, but of equal importance are the joints with which the actuators are attached to the base and the moving top plate. With size and rigidity in mind, our engineers came up with innovative spherical joints that are not only very simple, but also enhance the stiffness of the hexapod.



As a standard feature, the HXP50-ELEC controller allows the user to choose a virtual pivot point in space for all rotations.

The HXP50-ELEC controller accurately masters the synchronized transformations from Cartesian input coordinates to the motion of the Hexapod legs. In addition, the HXP50-ELEC provides advanced features including instrument grade I/O's, hardware based input triggers, event triggers, high-speed on-the-fly data acquisition, fast TCP/IP communication, and integrated TCL programming language for on-board processes. All these features improve accuracy and throughput, making the programmer's life much easier.

The main differentiator of the HXP50-MECA, as with the other Newport hexapods, is the ability to program two pivot points represented by the Tool and Work coordinate systems. The Tool CS moves with the top plate and the Work is stationary. Imagine a machine tool where one can adjust the orientation of both the cutting tool and workpiece or in photonics, the optical beam and the sample. Incremental displacements are possible in either one in user-friendly Cartesian coordinates, and positions can be easily switched from one system to the other by a function call or by numerical input in the HXP's web site interface. These powerful functions are a completely new way of mastering Hexapod motions without the need for complex external coordinate transformations.

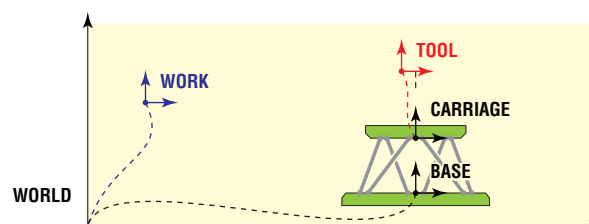

Absolute moves and positions are defined in the work coordinate system. Incremental moves can be done in the tool or in the work coordinate systems.

3.0 Characteristics

3.1 Mechanical Specifications

| | X | Y | Z | Θ_x | Θ_y | Θ_z |
|--|--------------------|--------------------|-----------------|------------------------|------------------------|------------------------|
| Travel range ⁽¹⁾ | ±17 mm | ±15 mm | ±7 mm | ±9° | ±8.5° | ±18° |
| MIM, Minimum incremental motion | 0.1 µm | 0.1 µm | 0.05 µm | 0.05 mdeg | 0.05 mdeg | 0.1 mdeg |
| Uni-directional repeatability, typical | 0.2 µm | 0.2 µm | 0.1 µm | 0.1 mdeg | 0.1 mdeg | 0.2 mdeg |
| Bi-directional repeatability, typical ⁽²⁾ | 1.2 µm (0.4 µm) | 1.2 µm (0.4 µm) | 0.6 (0.2 µm) | 0.6 mdeg (0.2 mdeg) | 0.6 mdeg (0.2 mdeg) | 1.2 mdeg (0.4 mdeg) |
| Max. speed | 14 mm/s | 12 mm/s | 5 mm/s | 6 °/s | 6 °/s | 15 °/s |
| Stiffness ⁽³⁾ | 2 N/µm | 2 N/µm | 25 N/µm | — | — | — |

¹⁾ Travel ranges are interdependent. The listed values are max. travels per axis when all other axis are in their centered position (Height = 151 mm for Z).

²⁾ With standard compensation (with hysteresis compensation).

³⁾ Stiffness depends on Hexapod position. Values are given for all axis in their centered position.




CAUTION


To reach specifications stated, stages must be fixed on a plane surface with a flatness of 50 µm or better.

3.2 Max. Cantilever Distance of the Load


3.2.1 Horizontal Base Plate

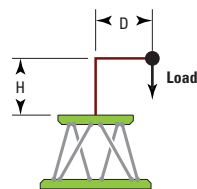


Load Position: D = 0 to 120 mm; H = 0 to 150 mm

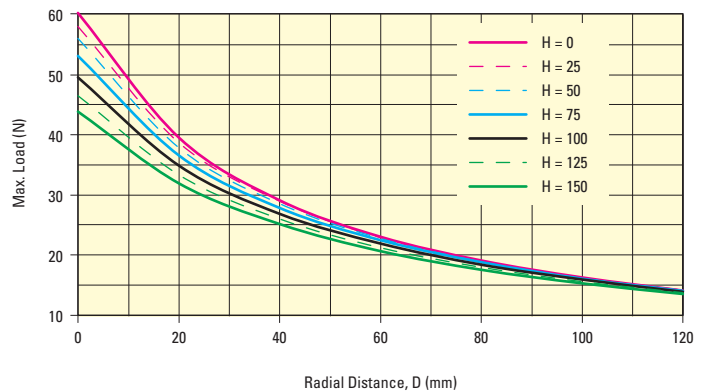


3.2.2 Base Plate Upside-Down

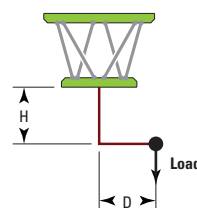


Load Position: D = 0 to 120 mm; H = 0 to 150 mm

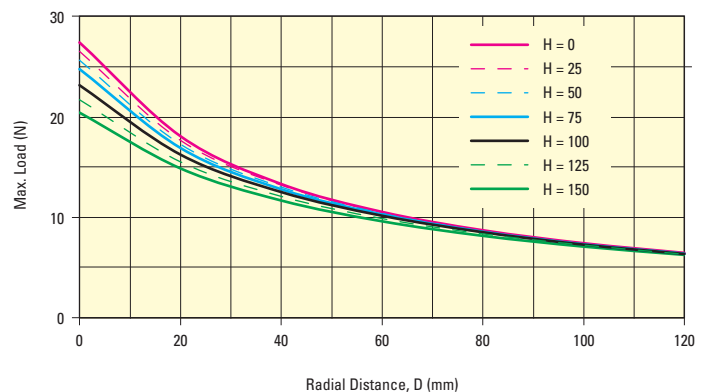


3.2.3 Vertical Base Plate

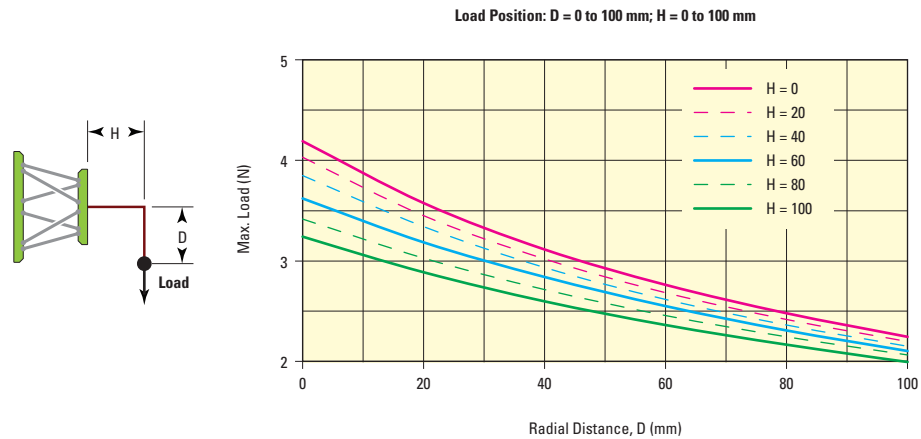

3.2.4 Base Plate at Any Position

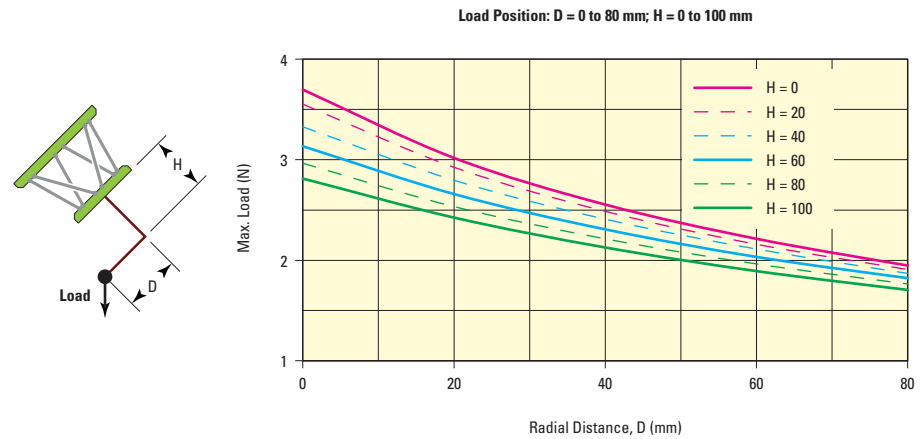

3.2.5 Horizontal Base Plate - Lateral Force

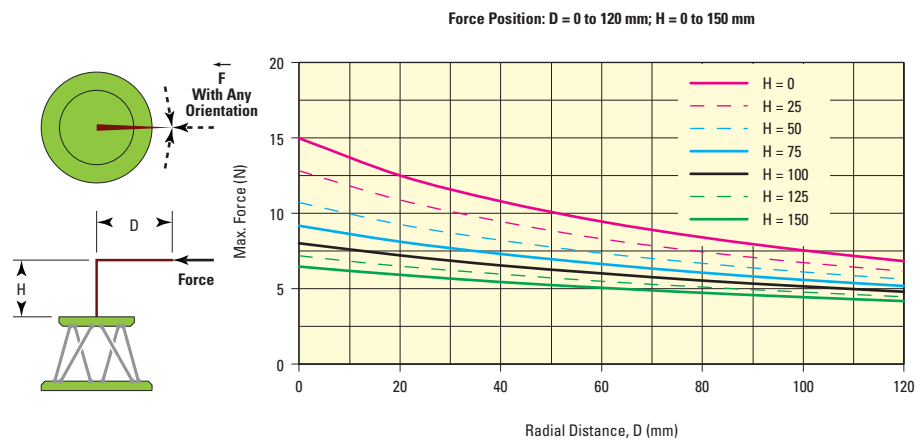

3.3 Weight

| | Weight [lb (kg)] |
|--------------------------|------------------|
| Hexapod (without cables) | 4.7 (2.15) |
| Cables | 2.8 (1.25) |

4.0 Connection to the HXP50-ELEC Controllers

4.1 Warnings on Controllers

Controllers are intended for use by qualified personnel who recognize shock hazards and are familiar with safety precautions required to avoid possible injury. Read the controller user's manual carefully before operating the instrument and pay attention to all written warnings and cautions.

WARNING

Disconnect the power plug under the following circumstances:

- If the power cord or any attached cables are frayed or damaged in any way.
- If the power plug is damaged in any way.
- If the unit is exposed to rain, excessive moisture, or liquids are spilled on the unit.
- If the unit has been dropped or the case is damaged.
- If you suspect service or repair is required.
- Whenever you clean the electronics unit.

CAUTION

To protect the unit from damage, be sure to:

- Keep all air vents free of dirt and dust.
- Keep all liquids away from the unit.
- Do not expose the unit to excessive moisture (85% humidity).
- Read this manual before using the unit for the first time.



WARNING

All attachment plug receptacles in the vicinity of this unit are to be of the grounding type and properly polarized.

Contact your electrician to check your receptacles.

WARNING

This product is equipped with a 3-wire grounding type plug.

Any interruption of the grounding connection can create an electric shock hazard.

If you are unable to insert the plug into your wall plug receptacle, contact your electrician to perform the necessary alterations to ensure that the green (green-yellow) wire is attached to earth ground.

WARNING

This product operates with voltages that can be lethal.

Pushing objects of any kind into cabinet slots or holes, or spilling any liquid on the product, may touch hazardous voltage points or short out parts.

4.2 Connection



WARNING

Always turn the controller's power OFF before connecting to a strut.

Anytime prior to power-up, each strut must be connected with the supplied cables to the corresponding 2 plugs located at the rear panel of the HXP50-ELEC controller. The SUB-D25 cable must be connected to the plug labeled “MOTOR”, the SUB-D15 cable to the plug labeled “ENCODER”.



**6 SUB-D25F CONNECTORS
FOR MOTOR CABLES**

**6 SUB-D15F CONNECTORS
FOR ENCODER CABLES**

WARNING

Each Hexapod strut label corresponds to the axis number of the HXP50-ELEC controller. This label is also indicated on connectors of the strut.



4.3 Cables

Each Hexapod is equipped with cables terminated with SUB-D25M connectors for motors and SUB-D15M connectors for encoders. These cables are directly connected to the HXP50-ELEC.

There are two cables per strut: 12 cables total to connect to the corresponding plugs on the controller.

CAUTION

Numbering of axes is important for both motor and encoder connections. A mismatch will cause fatal errors during initialization.

WARNING

These cables are shielded correctly. For a correct operation, make sure to lock connectors (ground continuity provided by the cable).

**WARNING**

Keep the motor cables at a safe distance from other electrical cables in your environment to avoid potential cross talk.

NOTE

On each connector, a sticker indicates the home preset and backlash on the SUB-D25M connector, and the resolution on the SUB-D15M connector.



These parameters have to be set into the HXP50-ELEC controller for each strut. This is done at our factory when the HXP50-MECA Hexapod is ordered with the HXP50-ELEC controller.




WARNING

The HXP50-MECA must be grounded via a threaded hole to avoid electrical disturbances generated by ground loops.



5.0

Dimensions

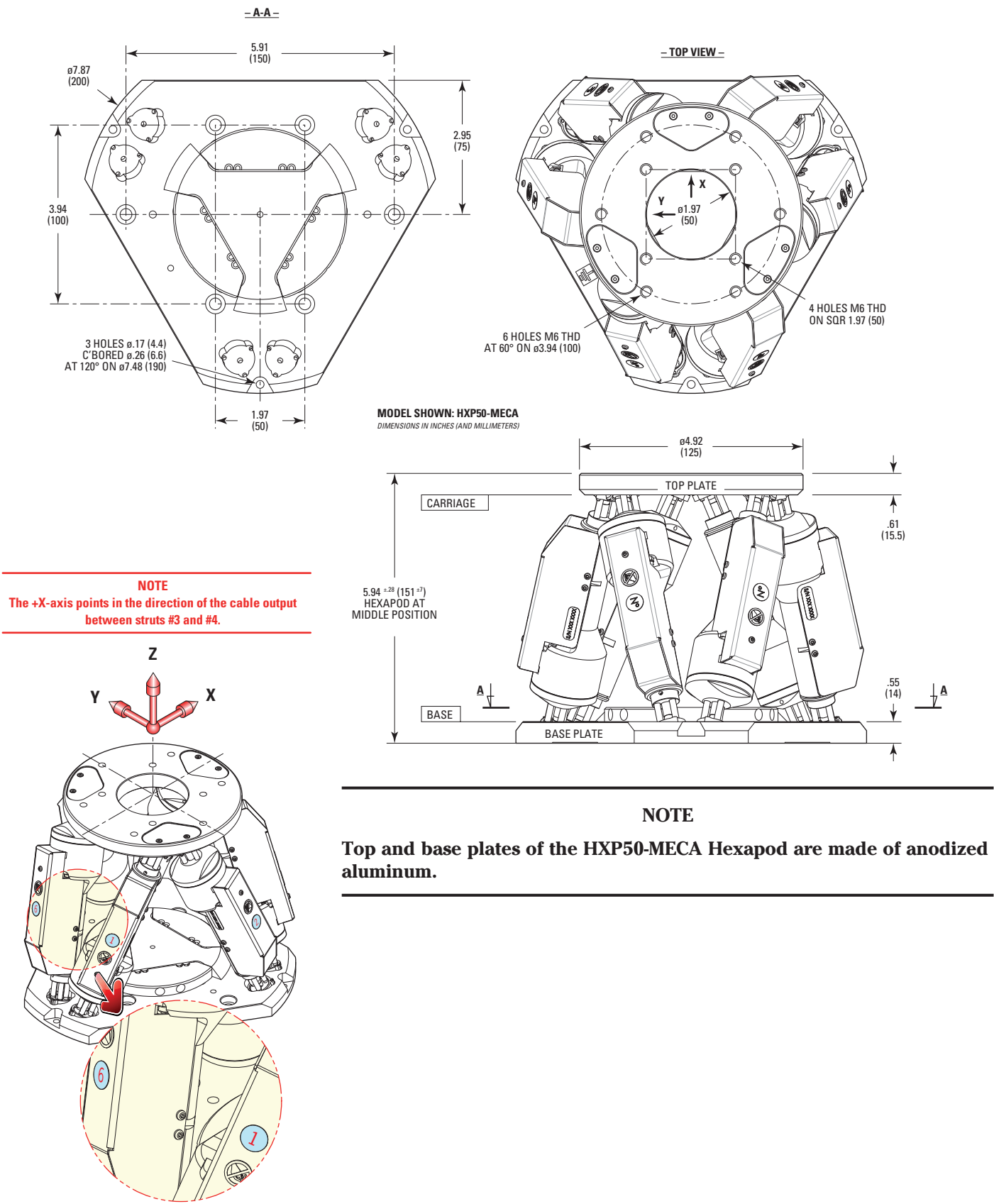

6.0 Maintenance

RECOMMENDATION

It is recommended to contact our After Sales Service which will know to define the appropriate maintenance for your application.

6.1 Maintenance

The Hexapod requires no particular maintenance. Nevertheless, this is a precision mechanical device that must be kept and operated with caution.

PRECAUTIONS

The Hexapod must be used or stocked in a clean environment, without dust, humidity, solvents or other substances.

RECOMMENDATION

It is recommended to return your Hexapod to Newport's After Sales Service after every 2000 hours of use for lubrication.

If your Hexapod is mounted on a workstation and cannot be easily removed, please contact Newport's After Sales Service for further instructions.

6.2 Repair

CAUTION



Never attempt to disassemble a component of the Hexapod that has not been covered in this manual.

To disassemble a non specified component can cause a malfunction of the stage.

If you observe a malfunction in your Hexapod, please contact us immediately to arrange for a repair.



CAUTION

Any attempt to disassemble or repair a Hexapod without prior authorization will void your warranty.

Service Form

Your Local Representative

Tel.: _____

Fax: _____

Name: _____

Return authorization #: _____

(Please obtain prior to return of item)

Company: _____

Address: _____

Date: _____

Country: _____

Phone Number: _____

P.O. Number:

Fax Number:

Item(s) Being Returned:

Model #: _____

Serial #: _____

Description: _____

Reasons of return of goods (please list any specific problems): _____



Newport®

Experience | Solutions

**Visit Newport Online at:
www.newport.com**

North America & Asia

Newport Corporation
1791 Deere Ave.
Irvine, CA 92606, USA

Sales

Tel.: (800) 222-6440
e-mail: sales@newport.com

Technical Support

Tel.: (800) 222-6440
e-mail: tech@newport.com

Service, RMAs & Returns

Tel.: (800) 222-6440
e-mail: service@newport.com

Europe

MICRO-CONTROLE Spectra-Physics S.A.S
9, rue du Bois Sauvage
91055 Évry CEDEX
France

Sales & Technical Support

Tel.: +33 (0)1.60.91.68.68
e-mail: france@newport.com

Service & Returns

Tel.: +33 (0)2.38.40.51.55